



180015144061



中国认可  
国际互认  
检测  
TESTING  
CNAS L2954

## Final Report

Report Number: SDWH-M202101057-4(E)

# Skin Irritation Test of Disposable Medical Protective Clothing

According to ISO 10993-10:2010  
0.9% Sodium Chloride Injection Extract

Sponsor: Wujiang Tutaike Textiles & Finishing Co.,Ltd

Address: No.1599,South 3rd Ring  
Road,Shengze,Wujiang,Suzhou,Jiangsu



**Sanitation & Environment Technology Institute, Soochow University**

Address: 199 Ren-Ai Road, Suzhou Industrial Park, Suzhou, Jiangsu 215123, P. R. China

Website: [www.sudatest.com](http://www.sudatest.com)

E-mail: [med@sudatest.com](mailto:med@sudatest.com)

Direct: +86 512 65880038

Free: 400 107 8828

## Content

<b>Supplementary Explanation</b> .....	<b>3</b>
<b>Quality Assurance Statement</b> .....	<b>4</b>
<b>GLP Compliance Statement</b> .....	<b>5</b>
<b>Verification Dates</b> .....	<b>5</b>
<b>Summary</b> .....	<b>6</b>
<b>Test Report</b> .....	<b>7</b>
<b>1 Purpose</b> .....	<b>7</b>
<b>2 Reference</b> .....	<b>7</b>
<b>3 Compliance</b> .....	<b>7</b>
<b>4 Identification of Test and Control Articles</b> .....	<b>7</b>
4.1 Test Article .....	7
4.2 Control Article.....	8
4.2.1 Negative Control .....	8
4.2.2 Positive Control.....	8
<b>5 Equipment and Reagents</b> .....	<b>8</b>
5.1 Equipment .....	8
5.2 Reagents .....	8
<b>6 Identification of Test System</b> .....	<b>8</b>
<b>7 Animal Care and Maintenance</b> .....	<b>9</b>
<b>8 Justification of Test System and Route of Administration</b> .....	<b>9</b>
<b>9 Experimental Design</b> .....	<b>9</b>
9.1 Preparation of Extracts.....	9
9.1.1 Pretreatment .....	9
9.1.2 Extraction .....	9
9.2 Experimental Procedure.....	9
9.3 Observation of Animals.....	10
9.4 Evaluation of Results .....	10
<b>10 Results</b> .....	<b>11</b>
<b>11 Conclusion</b> .....	<b>11</b>
<b>12 Record Storage</b> .....	<b>11</b>
<b>13 Confidentiality Agreement</b> .....	<b>11</b>
<b>14 Deviation Statement</b> .....	<b>11</b>
<b>Annex 1 Test Data</b> .....	<b>12</b>
<b>Annex 2 Photograph of Test Article</b> .....	<b>13</b>
<b>Annex 3 Information Provided by Sponsor</b> .....	<b>14</b>

## Supplementary Explanation

- (1) Please apply for rechecking within 15 days of receiving the report if there are any objections.
- (2) Any erasure or without special inspection and testing seal renders the report null and void.
- (3) The report is only valid when signed by the persons who edited, checked and approved it.
- (4) The results relate only to the articles tested.
- (5) The report shall not be reproduced except in full without the written approval of the institute.

## Quality Assurance Statement

The Quality Assurance Unit inspected/audited this study in compliance with the following GLP regulations:

Good Laboratory Practice (GLP) Regulation 21 CFR Part 58, U.S. Food and Drug Administration (FDA).

The Quality Assurance Unit conducted inspections on the following dates. The findings were reported to the Study Director and to the Testing Facility Management. The final report was reviewed by the Quality Assurance Unit. The final report accurately describes the test methods in accordance with standard operating procedures, and the results are consistent with raw data of non-clinical studies conducted according to the study protocol.

Inspections	Date of Inspection	Date Reported to Study Director	Date Reported to Testing Facility Management.
Study Protocol	2021-03-10	2021-03-10	2021-03-29
Study Procedure	2021-03-16	2021-03-16	2021-03-29
Raw Data	2021-03-29	2021-03-29	2021-03-29
Final Report	2021-03-29	2021-03-29	2021-03-29

Quality Assurance Unit: Zou Jing

Quality Assurance

2021-03-29

Date

## GLP Compliance Statement

This study was fully in accordance with the technical requirements of the study protocol.

This study was conducted in compliance with Good Laboratory Practice (GLP) Regulation 21 CFR Part 58, U.S. Food and Drug Administration (FDA).

### Verification Dates

Test Article Receipt	2021-03-08
Protocol Effective Date	2021-03-10
Technical Initiation Date	2021-03-11
Technical Completion Date	2021-03-19
Final Report Completion Date	2021-04-07

Edited by:

Xu Yueque2021-03-26

Date

Reviewed by:

Jimingwei2021-04-07

Date

Study Director

Approved by:

Fangyingyi2021-04-07

Date

Authorized Signatory

Sanitation &amp; Environment Technology Institute, Soochow University



# Summary

## 1 Test Article

<b>Test Article Name</b>	Disposable Medical Protective Clothing
<b>Manufacturer</b>	Wujiang Tutaike Textiles & Finishing Co.,Ltd
<b>Address</b>	No.1599,South 3rd Ring Road,Shengze,Wujiang,Suzhou,Jiangsu
<b>Model</b>	Not supplied by sponsor (N/S)
<b>Lot/Batch</b>	TTK-20200818

## 2 Main Reference

ISO 10993-10:2010 Biological evaluation of medical devices — Part 10: Tests for irritation and skin sensitization

## 3 Test Method

The extract of test article was evaluated for skin irritation. With ISO 10993-10:2010 Biological evaluation of medical devices — Part 10: Tests for irritation and skin sensitization.

Study protocol number: SDWH-PROTOCOL- GLP-M202101057-4.

## 4 Conclusion

The test result showed that the response of the test article extract was categorized as negligible under the test condition.

# Test Report

## 1 Purpose

The extract of test article was evaluated for skin irritation and extrapolating the results to humans, but it does not establish the actual risk of irritation.

## 2 Reference

ISO 10993-10:2010 Biological evaluation of medical devices — Part 10: Tests for irritation and skin sensitization

ISO 10993-12:2021 Biological evaluation of medical devices — Part 12: Sample preparation and reference materials

ISO 10993-2:2006 Biological evaluation of medical devices — Part 2: Animal welfare requirements

## 3 Compliance

Good Laboratory Practice Regulations, 21 CFR, Part 58.

ISO/IEC 17025:2017 General requirements for the competence of testing and calibration laboratories (CNAS—CL01 Accreditation criteria for the competence of testing and calibration laboratories) China National Accreditation Service for Conformity Assessment LABORATORY ACCREDITATION CERTIFICATE Registration No. CNAS L2954.

RB/T 214—2017 Competence assessment for inspection body and laboratory mandatory approval—General requirements for inspection body and laboratory Certification and Accreditation Administration of the People's Republic of China INSPECTION BODY AND LABORATORY MANDATORY APPROVAL Certificate No. CMA 180015144061.

## 4 Identification of Test and Control Articles

### 4.1 Test Article

Test Article Name	Disposable Medical Protective Clothing
Manufacturer	Wujiang Tutaike Textiles & Finishing Co.,Ltd
Address	No.1599,South 3rd Ring Road,Shengze,Wujiang,Suzhou,Jiangsu
Test Article Initial State	Sterile, EO
CAS Number	N/S
Model	N/S
Size	N/S
Lot/Batch	TTK-20200818
Raw Material	PP+PE
Packaging Material	Paper-plastic bag
Physical State	Solid
Color	White
Density	N/S
Stability	N/S
Solubility	N/S
Storage Condition	Room temperature
Intended Use	N/S
Additional Information	N/S

The information about the test article was supplied by the sponsor wherever applicable.

The Sponsor is responsible for all test article characterization data as specified in the GLP regulations.

## 4.2 Control Article

### 4.2.1 Negative Control

Name: 0.9% sodium chloride injection (SC)  
 Manufacturer: Guangxi Yuyuan Pharmaceutical Co., Ltd.  
 Size: 500mL  
 Lot/ Batch#: H20110801  
 Physical State: Liquid  
 Color: Colourless  
 Storage Condition: Room Temperature

### 4.2.2 Positive Control

Name: sodium dodecyl sulfate  
 Manufacturer: Ron reagent  
 Size: 500g  
 Lot/ Batch#: RH178474  
 Physical State: Powder  
 Color: White  
 Storage Condition: Room Temperature  
 Solvent: 0.9% sodium chloride injection (SC)  
 Concentration: 20%  
 Date prepared: 2020-12-29

## 5 Equipment and Reagents

### 5.1 Equipment

Equipment Name	Equipment Number	Calibration Expire
Electronic Scale	SDWH2436	2021-05-21
Horizontal Large Capacity Constant Temperature Vibrator	SDWH2718	2021-08-10
Steel straight scale	SDWH463	2021-07-06
Vertical pressure steam sterilizer	SDWH2097	2022-03-09

### 5.2 Reagents

Reagent Name	Manufacturer	LOT
0.9% sodium chloride injection (SC)	Guangxi Yuyuan Pharmaceutical Co., Ltd.	H20110801
Sodium dodecyl sulfate (SDS)	Ron reagent	RH178474

## 6 Identification of Test System

Species: New Zealand white Rabbit (single strain).  
 Number: 3  
 Sex: Female  
 Weigh: Initial body weight not less than 2kg  
 Health status: Healthy, not previously used in other experimental procedures, young adult, nulliparous and not pregnant.  
 Housing: Animals were housed in cages identified by a card indicating the lab number, test code and first treatment date.  
 Animal identification: Stain with dyeing liquid  
 Cages: Stainless steel cage



Acclimation Period: 7 days under the same conditions as for the actual test

## 7 Animal Care and Maintenance

Animal purchase: Provided by Suzhou Experimental Animal Sci-tech Co., Ltd. <Permit Code: SCXK (SU) 2020-0007>

Bedding: NA

Feed: Rabbit Diet, Suzhou Experimental Animal Sci-tech Co., Ltd.

Water: Drinking water met the Standards for Drinking Water Quality GB 5749-2006

Animal room temperature: 18-26°C

Animal room relative humidity: 30%-70%

Lights: 12 hours light/dark cycle, full-spectrum lighting

Personnel: Associates involved were appropriately qualified and trained.

Selection: Only healthy, previously unused animals were selected.

There were no known contaminants present in the feed, water expected to interfere with the test data.

## 8 Justification of Test System and Route of Administration

The rabbit is specified as an appropriate animal model for evaluating potential skin irritants by the current testing standards. Positive control sodium dodecyl sulfate has been substantiated at SDWH with this method. See table 3.

The patches (about 2.5cm×2.5cm) which moistened by test article extract, and directly applying to the rabbit skin is considered to be the best mean of contact.

## 9 Experimental Design

### 9.1 Preparation of Extracts

#### 9.1.1 Pretreatment

No pretreatment required.

#### 9.1.2 Extraction

Under aseptic conditions, samples were taken according to the sampling method (Randomly sampling of the main material part). Calculate the standard surface area of the contact part (single side) of each sample. The extraction was performed with agitation in closed inert containers according to the extraction ratio listed in the following table (sample: extraction vehicle). The extraction vehicle was SC.

Test Period	Actual Sampling	Extract Procedure			Final Extract
		Extract Ratio	Extraction volume	Condition	
polar test extract	Surface area 120 cm <sup>2</sup>	6 cm <sup>2</sup> : 1 mL	20.0 mL	50°C, 72 h	Clear
polar negative control	/	/	10.0 mL		Clear


The state of the extract did not change after extraction. The extract was stored at room temperature, and tested within 24 h, without the process of adjusting its pH value, filtering, centrifuging, diluting, etc. The vehicle (without the test article) was similarly prepared to serve as the control.

### 9.2 Experimental Procedure

Use the rabbits with healthy intact skin. Fur was generally clipped within 4-24 h of testing on the backs of the rabbits, a sufficient distance on both sides of the spine for application and observation of all test sites (approximately 10cm×15 cm).

Apply 0.5 mL extract (s) of test article or control to 2.5 cm×2.5 cm absorbent gauze patches, and then apply the patch soaked with the extract of test article or control directly to the skin on each

side of each rabbit as shown in Figure 1, and then wrap the application sites with a bandage (semi-occlusive or occlusive) for a minimum of 4 h. At the end of the contact time, remove the dressing and washing with lukewarm water or other suitable nonirritating solvent and careful drying.

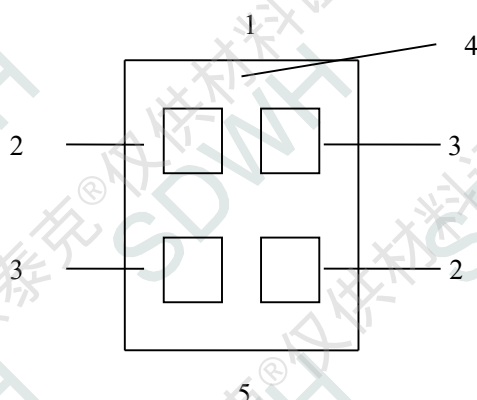

1- Cranial end, 2- Test site, 3- Control site, 4- Clipped dorsal region, 5- Caudal end

**Figure1 Location of skin application sites**

### 9.3 Observation of Animals

Describe and score the skin reaction for erythema and oedema according to the scoring system given in Table 1 for each application site at each time interval. Record the appearance of each application site at  $(1\pm 0.1)$  h,  $(24\pm 2)$  h,  $(48\pm 2)$  h and  $(72\pm 2)$  h following removal of the patches.

**Table 1 — Scoring system for skin reaction**

Reaction	Irritation score
<b>Erythema and Eschar Formation</b>	
No erythema	0
Very slight erythema (barely perceptible)	1
Well-defined erythema	2
Moderate erythema	3
Severe erythema (beet redness) to eschar formation preventing grading of erythema	4
<b>Oedema Formation</b>	
No edema	0
Very slight edema (barely perceptible)	1
Well-defined edema (edges of area well-defined by definite raising)	2
Moderate edema (raised approximately 1mm)	3
Severe edema (raised more than 1mm and extending beyond exposure area)	4
Maximal possible score for irritation	8
Other adverse changes at the skin sites shall be recorded and reported.	

### 9.4 Evaluation of Results

Use only  $(24\pm 2)$  h,  $(48\pm 2)$  h and  $(72\pm 2)$  h observations for calculation.

After the 72 h grading, all erythema grades plus oedema grades  $(24\pm 2)$  h,  $(48\pm 2)$  h and  $(72\pm 2)$  h were totalled separately for each test sample and blank for each animal. The primary irritation score for an animal was calculated by dividing the sum of all the scores by 6 (two test/observation sites, three time points).

To obtain the primary irritation index for the test article, add all the primary irritation scores of the individual animals and divide by the number of animals.

When blank or negative control is used, calculate the primary irritation score for the controls and

subtract that score from the score using the test material to obtain the primary irritation score. The primary irritation index (PII) for the test article was evaluated according to Table 2.

**Table 2 — Primary or cumulative irritation index categories in a rabbit**

Mean score	Response category
0~0.4	Negligible
0.5~1.9	Slight
2~4.9	Moderate
5~8	Severe

## 10 Results

All animals were survived and no abnormal signs were observed during the study. According to what observed, the response of skin on testing side did not exceed that on the control side. Thus, the primary irritation index for the test article was calculated to be 0. See table 4.

## 11 Conclusion

The test result showed that the response of the test article extract was categorized as negligible under the test condition.

## 12 Record Storage

All raw data pertaining to this study and a copy of the final report are to be retained in designated SDWH archive.

## 13 Confidentiality Agreement

Statements of confidentiality were as agreed upon prior to study initiation.

## 14 Deviation Statement

There were no deviations from the approved study protocol which were judged to have any impact on the validity of the data.

## Annex 1 Test Data

**Table 3** Positive control

Extract	Rabbit No.	Group	Reaction	Interval (hours): score=left site/right site		
				24±2h	48±2h	72±2h
SC	1	Positive Control	Erythema	2/2	3/3	3/3
			Oedema	3/3	4/3	4/4
		Negative Control	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
SC	2	Positive Control	Erythema	2/2	3/2	3/3
			Oedema	3/3	4/4	4/4
		Negative Control	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
SC	3	Positive Control	Erythema	2/2	2/3	3/3
			Oedema	3/3	4/3	4/4
		Negative Control	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
The primary irritation score.					6.1	

Note: Positive control performed once every six months, see SDWH-M202007070-1(Completed Date: 2021-01-01).

**Table 4** Test Results of Dermal Observations

Extract	Rabbit No.	Group	Reaction	Interval (hours): score=left site/right site		
				24±2h	48±2h	72±2h
SC	1	Test Article	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
		Negative Control	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
SC	2	Test Article	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
		Negative Control	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
SC	3	Test Article	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
		Negative Control	Erythema	0/0	0/0	0/0
			Oedema	0/0	0/0	0/0
The primary irritation score.					0	

## Annex 2 Photograph of Test Article



## **Annex 3 Information Provided by Sponsor**

### **1 Production Process**

Not supplied by sponsor.

### **2 Other Information**

Not supplied by sponsor.

---

End of Report